

(12) INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(19) World Intellectual Property Organization
International Bureau



(43) International Publication Date
9 January 2003 (09.01.2003)

PCT

(10) International Publication Number
WO 03/003792 A1


- (51) International Patent Classification⁷: **H04R 25/00**
- (21) International Application Number: PCT/DK02/00454
- (22) International Filing Date: 28 June 2002 (28.06.2002)
- (25) Filing Language: English
- (26) Publication Language: English
- (30) Priority Data:
PA 2001 01016 28 June 2001 (28.06.2001) DK
- (71) Applicant (for all designated States except US): **OTICON A/S [DK/DK]**; Strandvejen 58, DK-2900 Hellerup (DK).
- (72) Inventors; and
- (75) Inventors/Applicants (for US only): **LUNDH, Peter [DK/DK]**; Oticon A/S, Strandvejen 58, DK-2900 Hellerup (DK). **SCHUM, Don [US/US]**; Oticon A/S, Strandvejen 58, DK-2900 Hellerup (DK).
- (74) Agent: **CHRISTENSEN, Mikael, T.**; c/o Oticon A/S, Strandvejen 58, DK-2900 Hellerup (DK).
- (81) Designated States (national): AE, AG, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, BZ, CA, CH, CN, CO, CR, CU, CZ, DE, DK, DM, DZ, EC, EE, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, MZ, NO, NZ, OM, PH, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TN, TR, TT, TZ, UA, UG, US, UZ, VN, YU, ZA, ZM, ZW.
- (84) Designated States (regional): ARIPO patent (GH, GU, KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZM, ZW), Eurasian patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE, TR), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, GQ, GW, ML, MN, NE, SN, TD, TG).

Published:

- with international search report
- before the expiration of the time limit for amending the claims and to be republished in the event of receipt of amendments

[Continued on next page]

(54) Title: **HEARING AID FITTING**

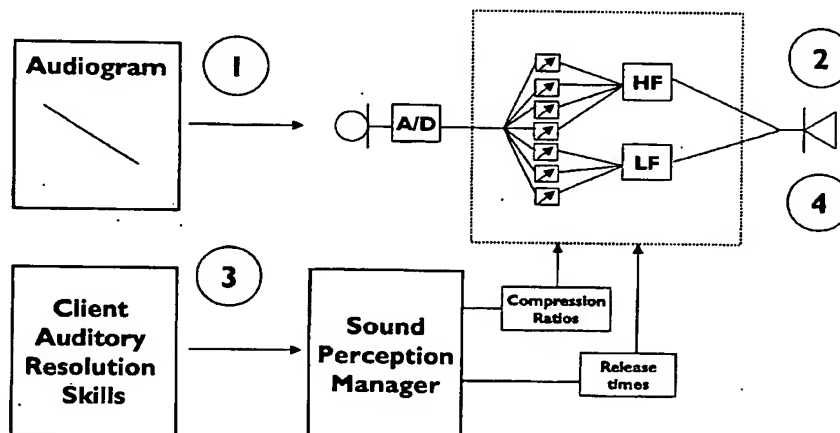

(57) Abstract: The invention relates to a fitting equipment for fitting a hearing aid, the equipment comprising a computer having a display, where the computer is provided with a fitting software adapted for controlling parameters of the hearing aid upon control of indicators in the software, where the indicators are visible on the display, where a default fitting rationale may be chosen and where in the fitting software at least two variant settings which are offset in relation to the default setting are provided and may be activated from the software by means of the computer, where in the first variant setting the release time in both the LF and HF channels are lengthened and the amount of gain for soft inputs is reduced in both the LF and HF channels and where in the second variant setting the release time in the HF is shortened and LF gain is increased in both the LF and HF channels.

WO 03/003792 A1

BEST AVAILABLE COPY

WO 03/003792 A1



For two-letter codes and other abbreviations, refer to the "Guidance Notes on Codes and Abbreviations" appearing at the beginning of each regular issue of the PCT Gazette.

TITLE

Hearing aid fitting

5

AREA OF THE INVENTION

The invention relates to the area of hearing aid fitting. More particularly the invention relates to the area of equipment for fitting hearing aid to the specific needs of hearing impaired individuals.

10

BACKGROUND OF THE INVENTION

Most modern hearing aids are programmable for adaptation to the hearing aid users needs. Rationales have been developed, which provides a good first approach to the fitting of the hearing aid to the user. The rationales are data sets specifying the transfer function or the gain of the hearing aid over a relevant frequency area

Some users with severe to profound hearing losses have however a special need for adaptation, which goes beyond what can be achieved by means of the traditionally developed rationales. These users often have a relatively normal upper hearing threshold, (normally designated the uncomfortable level (UCL), whereas the lower hearing threshold (HTL) has been significantly offset compared to that of a normal hearing person.

Speech in everyday environments can cover a 70 dB or greater dynamic range. From the softest elements of soft speech to the most intense elements of loud speech, the listener with normal hearing uses most of his/her dynamic range throughout the day. Multi-channel non-linear processing is designed to make use of most - if not all - of this dynamic range available to the patient with sensorineural hearing loss. However, as the hearing loss moves into the severe and then profound ranges, this core assumption needs to be modified, especially for some clients.

The assumption common to the fitting of other hearing losses up to this point is that the ear should be able to make reasonable use of amplified speech information, even when compressed. Given the dramatic nature of the damage in profound hearing loss, this assumption may not always be true. The damage pattern in some ears with profound
5 hearing loss may need a signal that is more linear whereas others may be able to make full use of a signal that is highly compressed.

The objective of the present invention is to provide equipment that can provide an improved initial setting of the hearing aid in order to achieve a faster and better fitting of
10 the hearing aid.

SUMMARY OF THE INVENTION

According to the invention the objective of the invention is achieved by means of the
15 fitting equipment as defined in claim 1.

According to the invention the objective is further achieved by means of a software for use in connection with a fitting equipment as defined in claim 2.

20

BRIEF DESCRIPTION OF THE DRAWINGS

FIG. 1 is a schematic view of a fitting software display depicting the invention.

25

DESCRIPTION OF A PREFERRED EMBODIMENT

Speech in everyday environments can cover a 70 dB or greater dynamic range. From the softest elements of soft speech to the most intense elements of loud speech, the listener with normal hearing uses most of his/her dynamic range throughout the day. Multi-
30 channel non-linear processing is designed to make use of most - if not all - of this dynamic range available to the patient with sensorineural hearing loss. However, as the hearing loss moves into the severe and then profound ranges, this core assumption needs to be modified, especially for some clients.

The assumption common to the fitting of other hearing losses up to this point is that the ear should be able to make reasonable use of amplified speech information, even when compressed. Given the dramatic nature of the damage in profound hearing loss, this assumption may not always be true. The damage pattern in some ears with profound hearing loss may need a signal that is more linear (called Type 1 clients) whereas others may be able to make full use of a signal that is highly compressed (called Type 3 clients).

10 Those ears that are sensitive to too much compression (called Type 1) likely have extensive damage throughout the inner ear and a minimal amount of remaining, intact inner hair cells. Signal resolution beyond the ability to detect pure-tones is very poor. These listeners probably can make use only of the peaks of the speech signal.

15 On the other hand, there are other ears with loss in the profound range that can actually make use of a signal that is even more compressed (called Type 3 clients). These ears probably have extensive outer hair cell loss and some amount of inner hair cell loss (accounting for the profound thresholds), but still have enough remaining signal resolution capacity to extract information from a compressed signal.

20

New Control with Three Settings

Based on these fundamental sound processing differences between clients, the fitting equipment or the software for the fitting equipment includes a tool in the trimmer panel.

25 This control with three alternative settings changes the manner in which the gain, amount of compression, and compression type are implemented in the fitting:

Mode 2: For clients called Type 2 with average auditory resolution. Default for fittings and consistent with general ASA2 principles, combining fast acting syllabic compression in the LF channel with slow acting Adaptive Gain in the HF channel. The goal is to provide a modestly compressed signal to ensure adequate speech audibility for moderate and loud speech. Softer speech signals can be made more audible with use of the Manual Override (volume control wheel).

Mode 1: For clients with low auditory resolution (called Type 1). Compared to Mode 2, the release time in both the LF and HF channels are lengthened, providing a greater linear effect for short term changes in the speech signal. In addition, the amount of gain for soft inputs is reduced in both the LF and HF channels by typically around 10 dB. This effect reduces the compression ratio. The total effect of all of these changes is to make the processing act more linear on a moment-to-moment basis, and to present primarily the speech peaks to the listener, but to continue to compensate for long-term changes in the overall input speech level.

Mode 3: For clients with high auditory resolution (called Type 3). Compared to Mode 2, the release time in the HF is shortened and the LF gain is increased in both the LF and HF channels (thus increasing the compression ratios). The sum total of these changes is to pack more of the speech signal within the remaining dynamic range, assuming that the ear has sufficient remain integrity to resolve this densely compacted signal.

By default, clients will be assumed to be Type 2, with fast acting compression applied in the low frequencies and slow acting compression applied in the highs.

From FIG. 1 a display appears schematically showing the activation buttons by which the client type may be selected.

The following tables depict a preferred embodiment of the settings, which are selected by selecting the client type

	SOFT GAIN	
Mode 1	Mode 2	Mode 3
-8 dB	X	+3 dB

	RELEASE TIMES	[ms]	
	Mode 1	Mode 2	Mode 3
LF	160	80	80
HF	640	320	80

From FIG. 1 the general principles appear. The hearing aid is initially set according to the actual hearing loss. This is Mode 2. Based on the hearing impaired individuals
5 auditory resolution skills the hearing aid professional is able to select one of the two possibilities of variant settings, which are Mode 1 or Mode 3. Hereby the compression ratios and the release times may be changed in a reliable manner without any discomfort for the hearing impaired.

CLAIMS

1. A fitting equipment for fitting a hearing aid, the equipment comprising a computer having a display, where the computer is provided with a fitting software adapted for controlling parameters of the hearing aid upon control of indicators in the software, where the indicators are visible on the display, where a default fitting rationale may be chosen and where in the fitting software at least two variant settings which are offset in relation to the default setting are provided and may be activated from the software by means of the computer, where in the first variant setting the release time in both the LF and HF channels are lengthened and the amount of gain for soft inputs is reduced in both the LF and HF channels and where in the second variant setting the release time in the HF is shortened and the LF gain is increased in both the LF and HF channels.
2. A fitting equipment according to claim 1, where the release times are increased by more than 50% in a first variant setting in both the LF and the HF area and where the release times are reduced by more than 50% in the HF area in a second variant setting.
3. A fitting equipment according to claim 1 or 2, where the soft gain is reduced by 3 dB, preferably around 8 dB in a first variant setting and where the soft gain is increased by more than 1 dB, preferably around 3 dB in a second variant setting.
4. A software for use in a fitting equipment for fitting a hearing aid, the equipment comprising a computer having a display, where the software is adapted for controlling parameters of the hearing aid upon control of indicators in the software, where the indicators are visible on the display, where a default fitting rationale may be chosen and where in the fitting software at least two variant settings which are offset in relation to the default setting are provided and may be activated from the software by means of the computer, where in the first variant setting the release time in both the LF and HF channels are lengthened and the amount of gain for soft inputs is reduced in both the LF and HF channels and

where in the second variant setting the release time in the HF is shortened and the LF gain is increased in both the LF and HF channels

- 5 5. A software according to claim 4, where the release times are increased by more than 50% in a first variant setting in both the LF and the HF area and where the release times are reduced by more than 50% in the HF area in a second variant setting.
- 10 6. A software according to claim 4 or 5, where the soft gain is reduced by 3 dB, preferably around 8 dB in a first variant setting and where the soft gain is increased by more than 1 dB, preferably around 3 dB in a second variant setting.

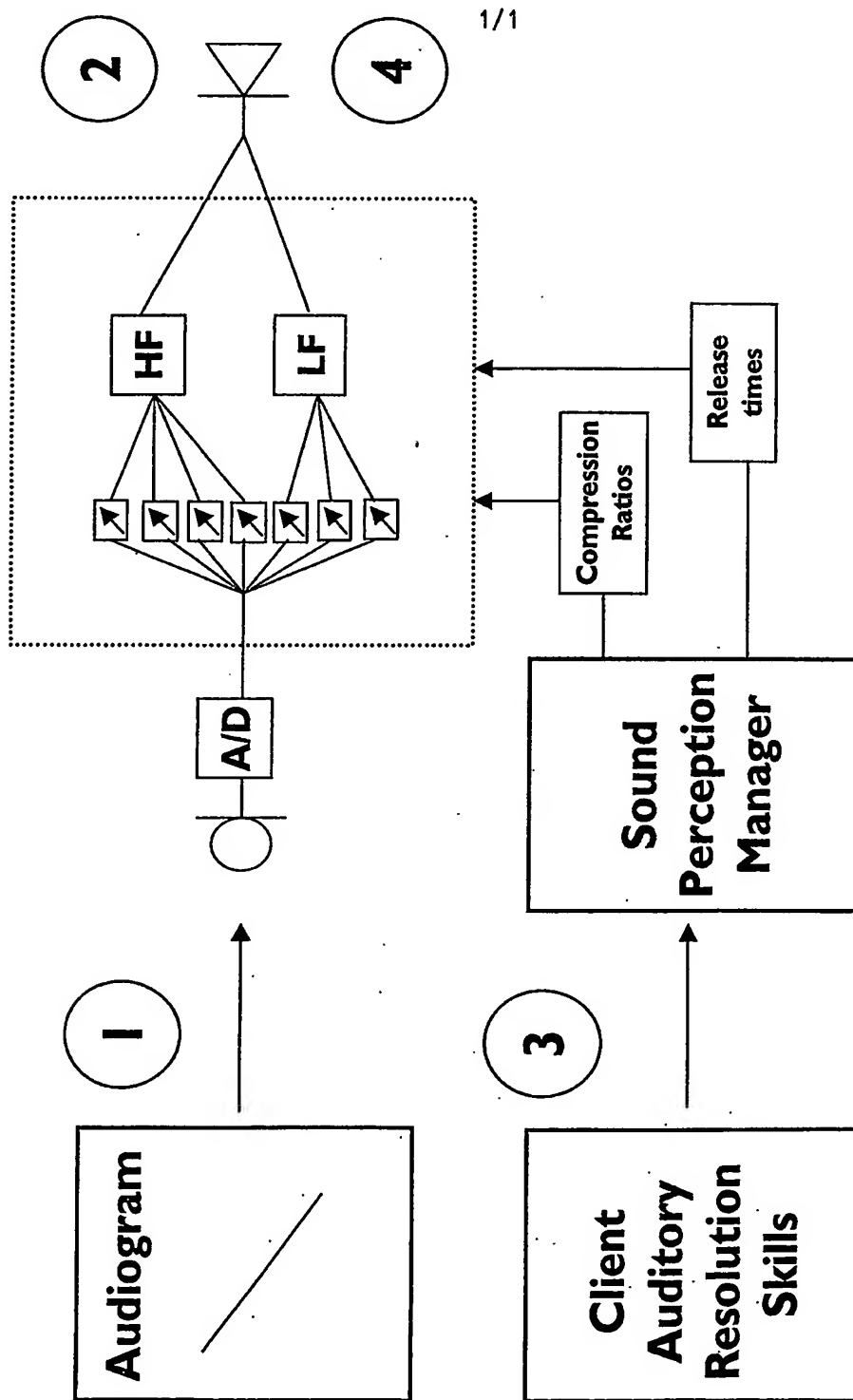


FIG. 1

INTERNATIONAL SEARCH REPORT

International Application No.

PCT/DK 02/00454

A. CLASSIFICATION OF SUBJECT MATTER

IPC 7 H04R25/00

According to International Patent Classification (IPC) or to both national classification and IPC

B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)

IPC 7 H04R

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practical, search terms used)

EPO-Internal, WPI Data

C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
A	US 6 201 875 B1 (DAVIS KEITH L. ET AL.) 13 March 2001 (2001-03-13) the whole document	1-6
A	US 5 991 417 A (TOPHOLM JAN) 23 November 1999 (1999-11-23) the whole document	1-6
A	US 5 838 801 A (ISHIGE RYUICHI ET AL) 17 November 1998 (1998-11-17) the whole document	1-6
P, A	WO 01 56331 A (TOEPHOLM & WESTERMANN ; LUDVIGSEN CARL (DK)) 2 August 2001 (2001-08-02) the whole document	1-6

-/-

☒ Further documents are listed in the continuation of box C.

☒ Patent family members are listed in annex.

* Special categories of cited documents:

"A" document defining the general state of the art which is not considered to be of particular relevance

"E" earlier document but published on or after the international filing date

"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)

"O" document referring to an oral disclosure, use, exhibition or other means

"P" document published prior to the international filing date but later than the priority date claimed

"T" later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention

"X" document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone

"Y" document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art

"&" document member of the same patent family

Date of the actual completion of the international search

22 October 2002

Date of mailing of the international search report

11.11.2002

Name and mailing address of the ISA

European Patent Office, P.B. 5818 Patentlaan 2
NL - 2280 HV Rijswijk
Tel. (+31-70) 340-2040, Tx. 31 651 epo nl,
Fax: (+31-70) 340-3016

Authorized officer:

HENRIK ANDERSSON /JA A

INTERNATIONAL SEARCH REPORT

International Application No.
PCT/DK 02/00454

C.(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
P, A	WO 02 39784 A (SHIM YOON JOO) 16 May 2002 (2002-05-16) the whole document	1-6

INTERNATIONAL SEARCH REPORT

Information on patent family members

International Application No

PCT/DK 02/00454

Patent document cited in search report		Publication date	Patent family member(s)	Publication date
US 6201875	B1	13-03-2001	WO 9948323 A2	23-09-1999
US 5991417	A	23-11-1999	WO 9635314 A1	07-11-1996
			AT 171833 T	15-10-1998
			AU 698105 B2	22-10-1998
			AU 2561695 A	21-11-1996
			DE 69505155 D1	05-11-1998
			DE 69505155 T2	15-04-1999
			DK 824845 T3	21-06-1999
			EP 0824845 A1	25-02-1998
			JP 11505077 T	11-05-1999
US 5838801	A	17-11-1998	JP 2904272 B2	14-06-1999
			JP 10174195 A	26-06-1998
			DK 142997 A	30-07-1998
WO 0156331	A	02-08-2001	AU 2831101 A	07-08-2001
			WO 0156331 A1	02-08-2001
WO 0239784	A	16-05-2002	AU 2414502 A	21-05-2002
			WO 0239784 A1	16-05-2002

**This Page is Inserted by IFW Indexing and Scanning
Operations and is not part of the Official Record**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ **BLACK BORDERS**
- ☐ **IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**
- ☐ **FADED TEXT OR DRAWING**
- ☐ **BLURRED OR ILLEGIBLE TEXT OR DRAWING**
- ☐ **SKEWED/SLANTED IMAGES**
- ☐ **COLOR OR BLACK AND WHITE PHOTOGRAPHS**
- ☐ **GRAY SCALE DOCUMENTS**
- ☐ **LINES OR MARKS ON ORIGINAL DOCUMENT**
- ☐ **REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**
- ☐ **OTHER:** _____

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.